


~ Atomic spectra and qualitative
Spectral analysis ~

Rad Emission - Atomic


Summary

- the aim of this laboratory is:
 - to observe the discrete spectra of two gases and to compare with the continuous spectrum of the natural light;
 - to calibrate a spectroscope;
 - to use the calibration curve for the determination of mercury spectral lines wavelength
- The different amounts of energy that an atomic electron is allowed are called energy levels. When an electron has the smallest allowed amount of energy, it occupies the lowest energy level. This level is called ground state.


Absorption

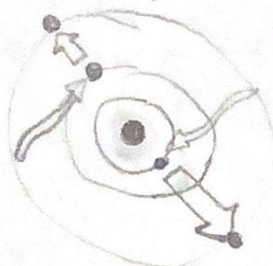

Emission

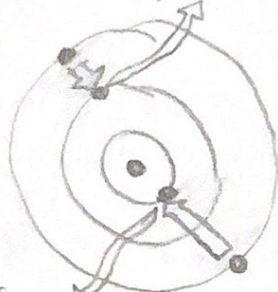

→ According to Bohr, the energy of an orbiting electron is an atom is the kinetic energy of the electron and the potential energy resulting from the Coulombian attraction force between the electron and nucleus.

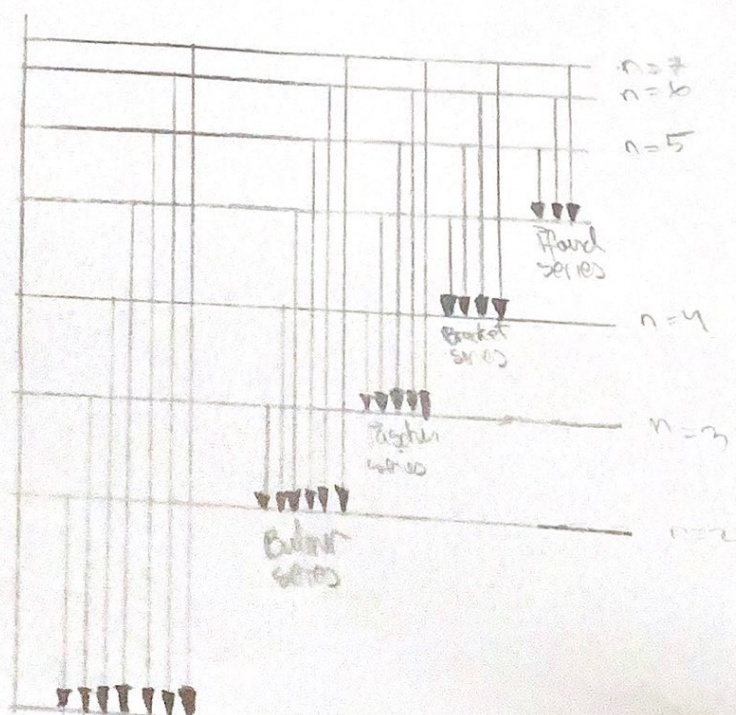

→ Einstein demonstrates that the light photon has a certain energy. Bohr postulated that a photon can be absorbed only if the energy of the photon is equal with the difference in the energy of the atomic electron levels.

$$h\nu = E_{\text{excited}} - E_{\text{ground}} \quad (1)$$

→ When the electron makes the transition to the ground state, a photon is emitted (Fig 1). The energy of the photon is equal to the energy difference the excited and the ground states. (Lyman series Fig 2)

→ The puzzle of the electron's position and motion around the nucleus in atom was clarified by studying the light emitted by atoms. The set of wave lengths of light emitted by an atom is called emission spectrum of that atom. When a body is heated it becomes incandescent.

Data : Neon							
Wave length $\lambda [Å]$	6402	6143	5945	5852	5400	5330	5330
line position [div]	660	590	595	580	575	540	532
line color	bright red	red orange	orange	yellow	bright green	soft green	light green

Data : Mercury		
color	line position	Wave length
Violet	405	4100
Blue	432	4400
Green	542	5600
Dark green	575	5800

$$E = 6,626 \cdot 10^{-34} \cdot \frac{3 \cdot 10^8}{4100 \cdot 10^{-10}} = 4,8 \cdot 10^{-19}$$

$$E = 4,5 \cdot 10^{-19} \text{ J}$$

$$E = 3,5 \cdot 10^{-19} \text{ J}$$

$$E = 3,3 \cdot 10^{-19} \text{ J}$$

$$\text{Violet: } 4,8 \cdot 10^{-19} : 1,6 \cdot 10^{-19} = 3 \text{ eV}$$

$$\text{Blue : } 4,5 \cdot 10^{-19} : 1,6 \cdot 10^{-19} = 2,8 \text{ eV}$$

$$\text{Green : } 3,5 \cdot 10^{-19} : 1,6 \cdot 10^{-19} = 2,18 \text{ eV}$$

$$\text{Dark Green: } 3,3 \cdot 10^{-19} : 1,6 \cdot 10^{-19} = 2,06 \text{ eV}$$

✓
fi

Lab 4

Atomic Spectra
Lab 3

$$\frac{E = h \cdot \nu}{\text{photon}} = \gamma \cdot \frac{1}{\lambda} = \gamma, \quad 1\text{eV} = 1.6 \cdot 10^{-19} \text{ J}$$

$$h = 6.626 \cdot 10^{-34} \text{ J} \cdot \text{s}, \quad c = 3 \cdot 10^8 \frac{\text{m}}{\text{s}}$$

Ne, Hg

$$\overset{\circ}{\text{A}} = 10^{-10} \text{ m}$$

$$\nu = \frac{c}{\lambda} = \frac{\overset{\text{speed of light}}{\text{m/s}}}{\text{m}} = \frac{1}{\lambda}$$

λ 700 nm 400 nm

FT-IR

UV

E =